

Controller template with DSA

For 2.14, 2.737

Instructor – Prof. David L. Trumper

Useful tips for LabVIEW:

- To open [LabVIEW in-context help](#) within LabVIEW, press Ctrl+H.
- To switch between [Front-panel and block diagram](#), press Ctrl+E.
- Always open the [top-level project file](#) with the “*.lvproj” extension to access the LV program.

Overview

- Controller template with up to 2 kHz loop rate.
- Dynamic Signal Analyzer to get the frequency response bode plots built in.
 - Plant
 - Loop return ratio
 - Etc.

Pinout Configuration

All connections will be made through the myRIO protoboard provided:

- AI 0 : Plant output being read by the myRIO
- AI 1 : Reference signal for the control loop
- AO 0 : Control effort signal to the plant


These connections are on the protoboard side of myRIO, called the Mini-system Port (MSP) or the “C” connector.

Front Panel

Loop timing

Set loop rate of each loop
(Restart program for
changes to take effect).

Also displays the
#Finished late counter
(refreshed after stopping
the program), which
shows the number of
times the particular loop
ran late. Useful to check
jitter in your discrete
control system.

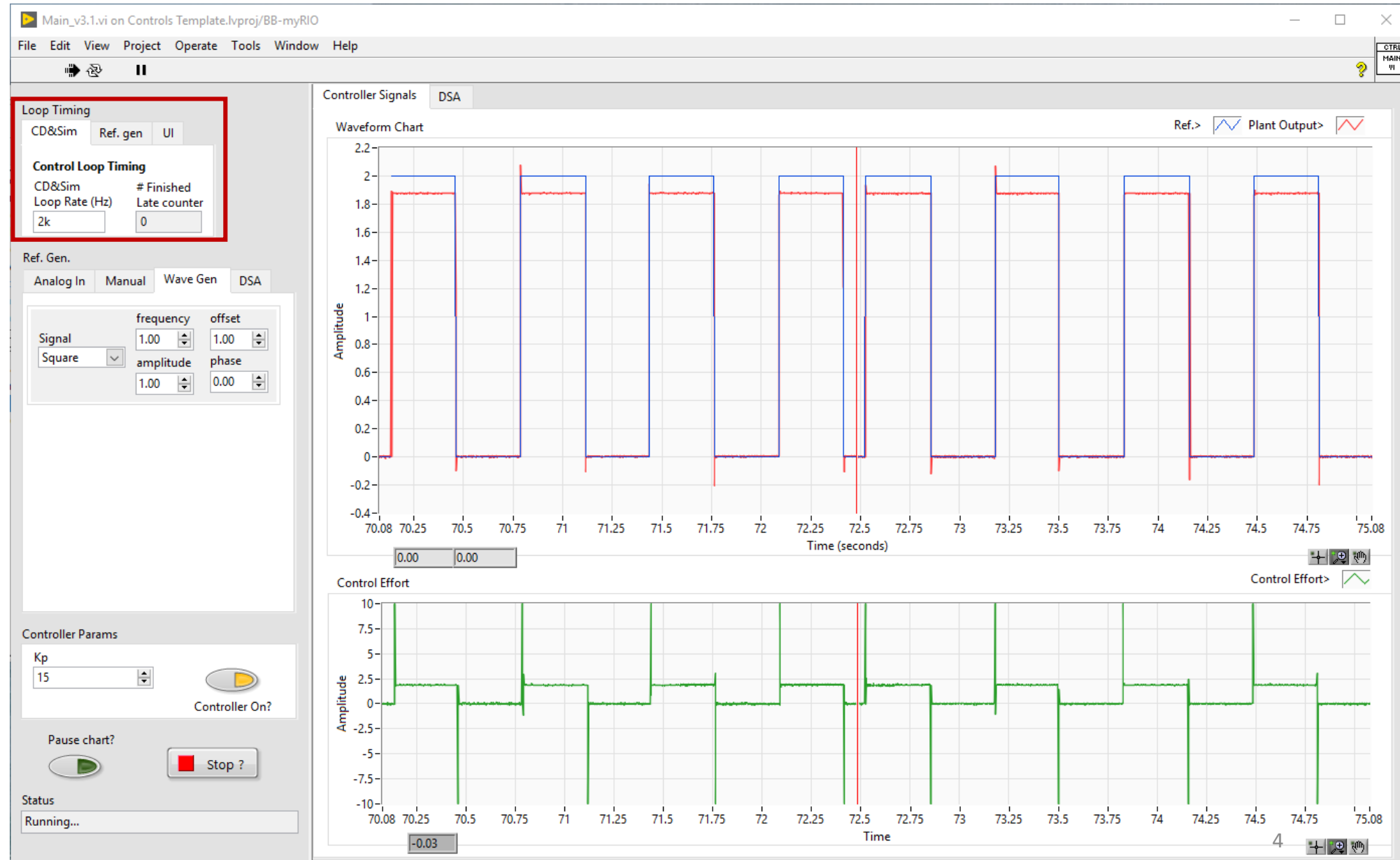


Front Panel


Ref. Gen

Various options of generating the reference signal.

- **Analog In:** Uses C/AI-1 as the input voltage signal.
- **Manual:** manually set DC value of the reference signal.
- **Wave Gen.:** use an in-built function generator to set the reference signal.
- **DSA:** Dynamic Signal Analyzer (used to get frequency responses).

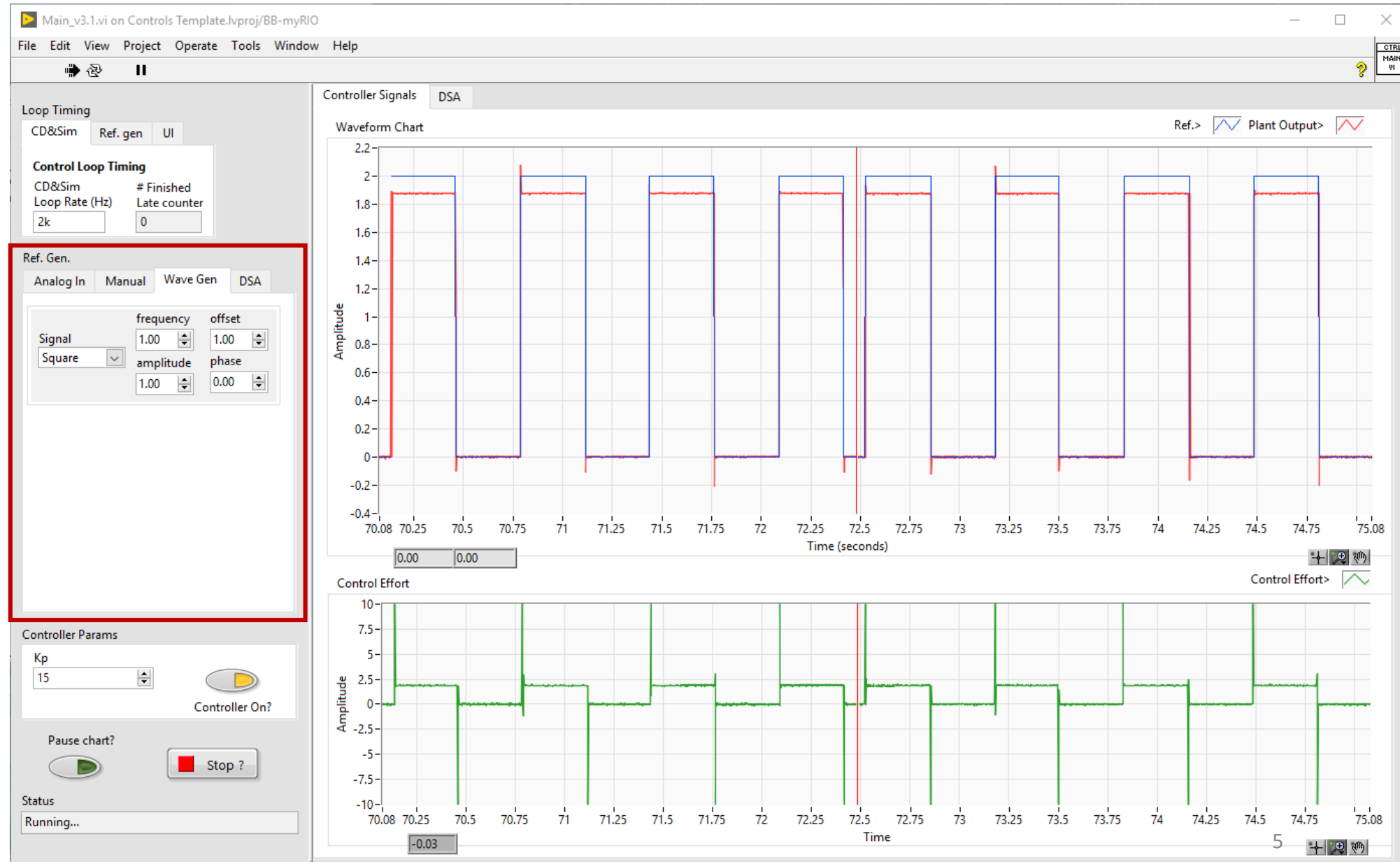

Front Panel

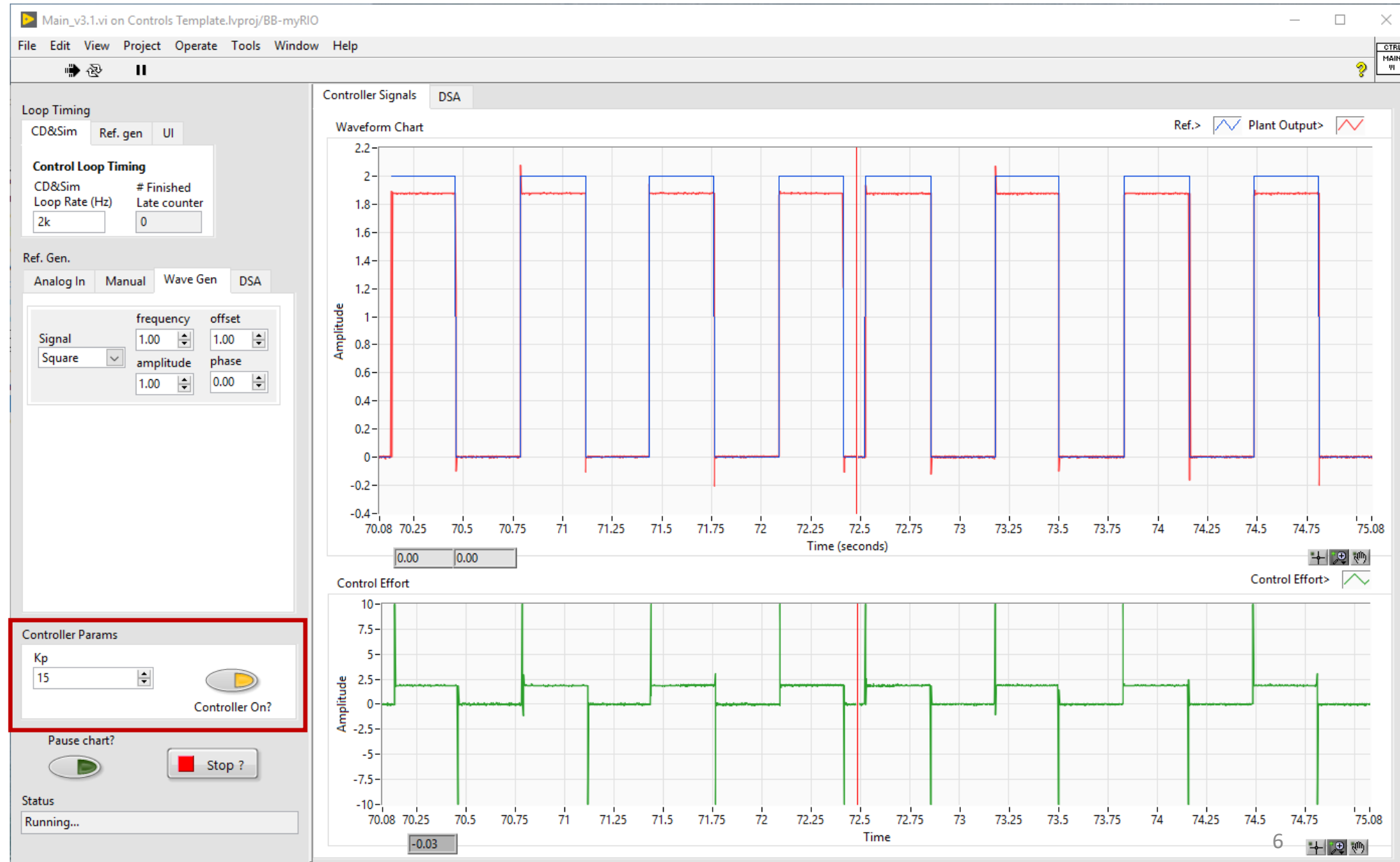


Controller Params

Set **Proportional Gain (Kp)** of the controller.

Controller On?:

If off, the controller output is forced to 0.

Front Panel

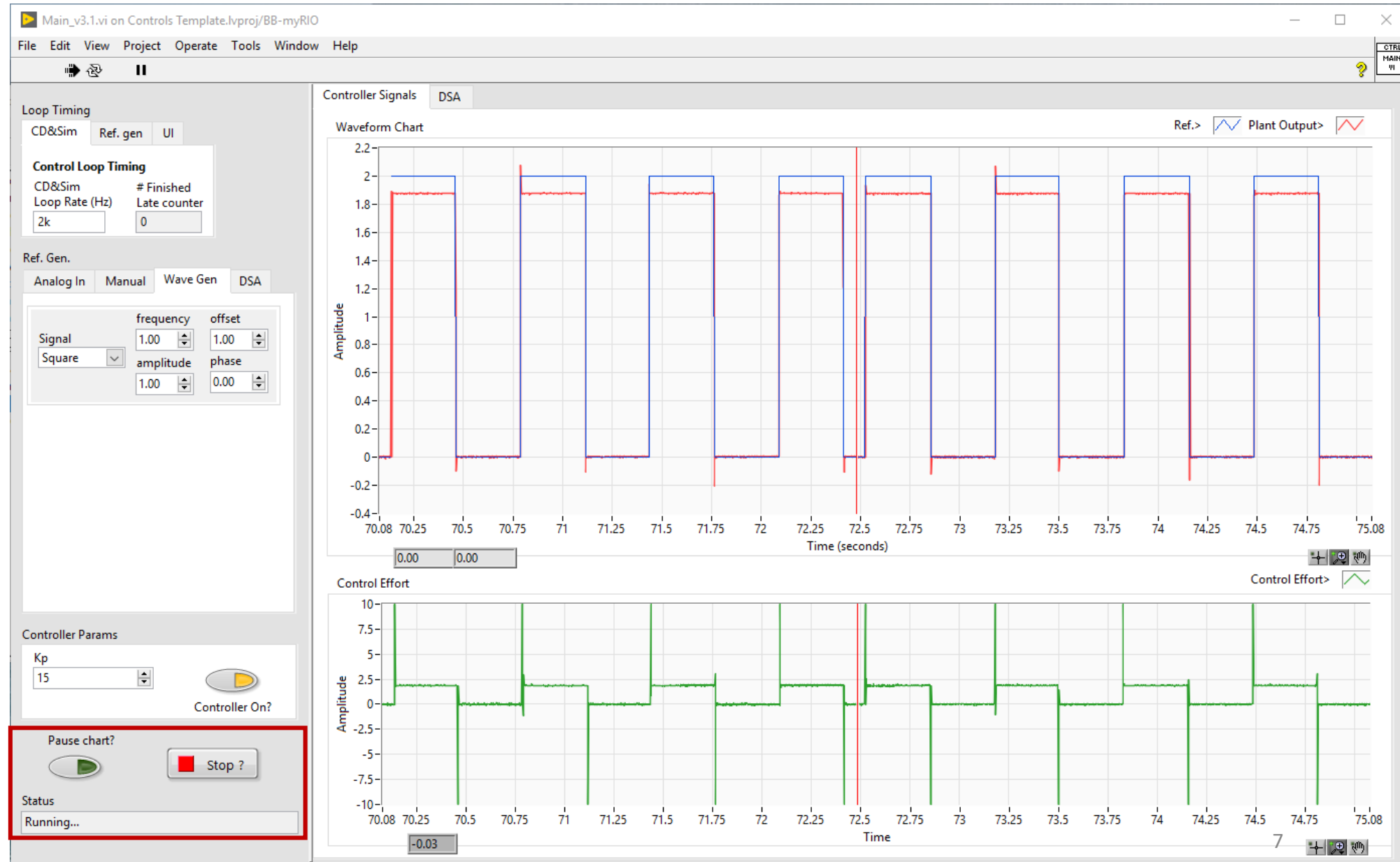

Program options

Pause chart:

Pauses the chart plotting to allow for zooming into the signals.

Stop?:

Stops the program

Status:

Displays current program status

Front Panel

Waveform tabs


Controller signals:

View reference, plant output and control-effort.

DSA:

View Dynamic Signal Analyzer (DSA) signals and frequency response measurement (after measurement is completed).


More on DSA later.

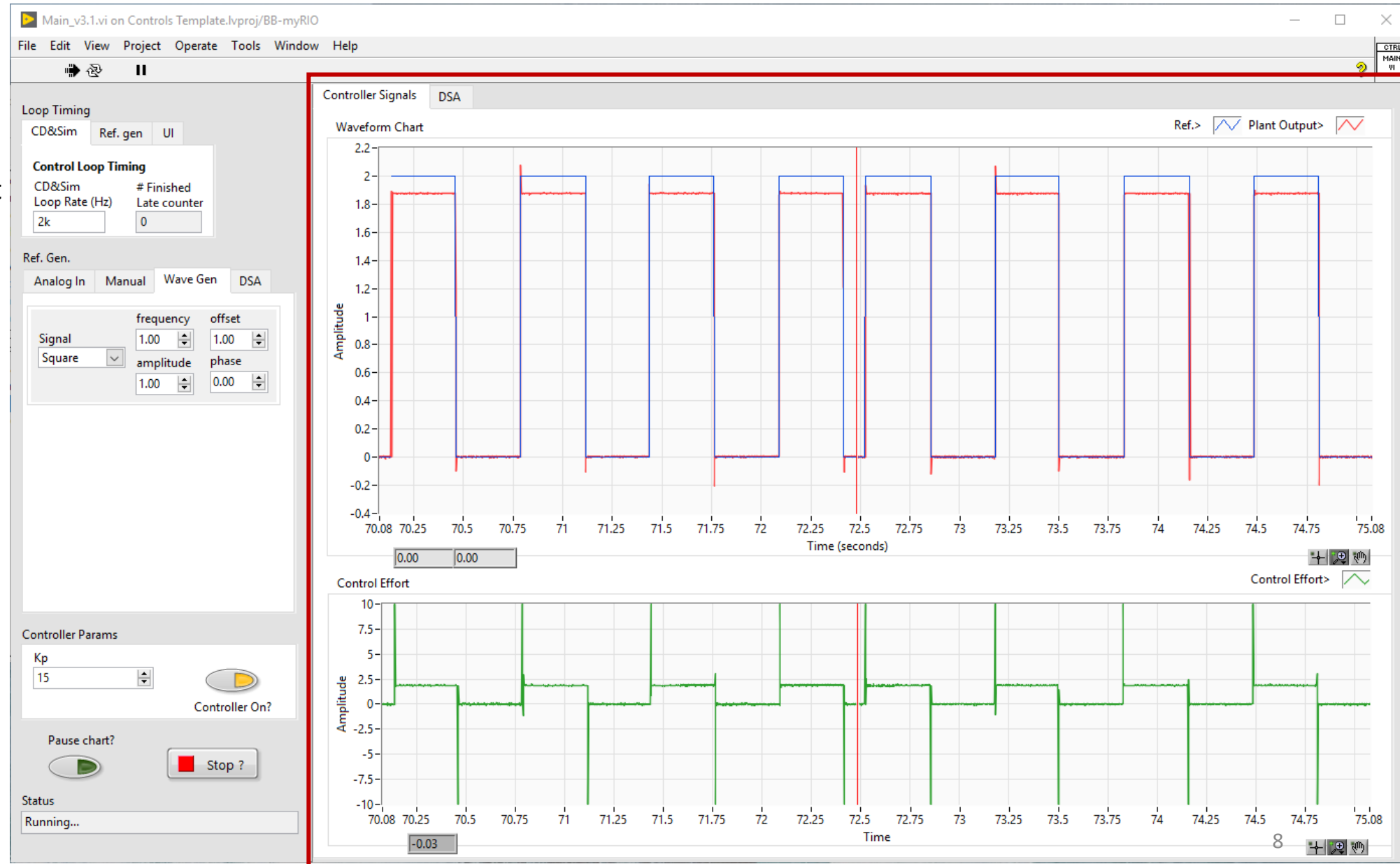


Block-diagram overview

- 3 parallel loops
 - User-Interface (UI) loop (non-deterministic)
 - [Control and Simulation Loop \(CD&Sim\)](#) – deterministic
 - Reference Generation / DSA loop – deterministic
- Inter-loop communication using [local variables](#) and [RT-FIFO](#) variables.
- [State machine architecture](#) for User-Interface.

Block diagram - controller

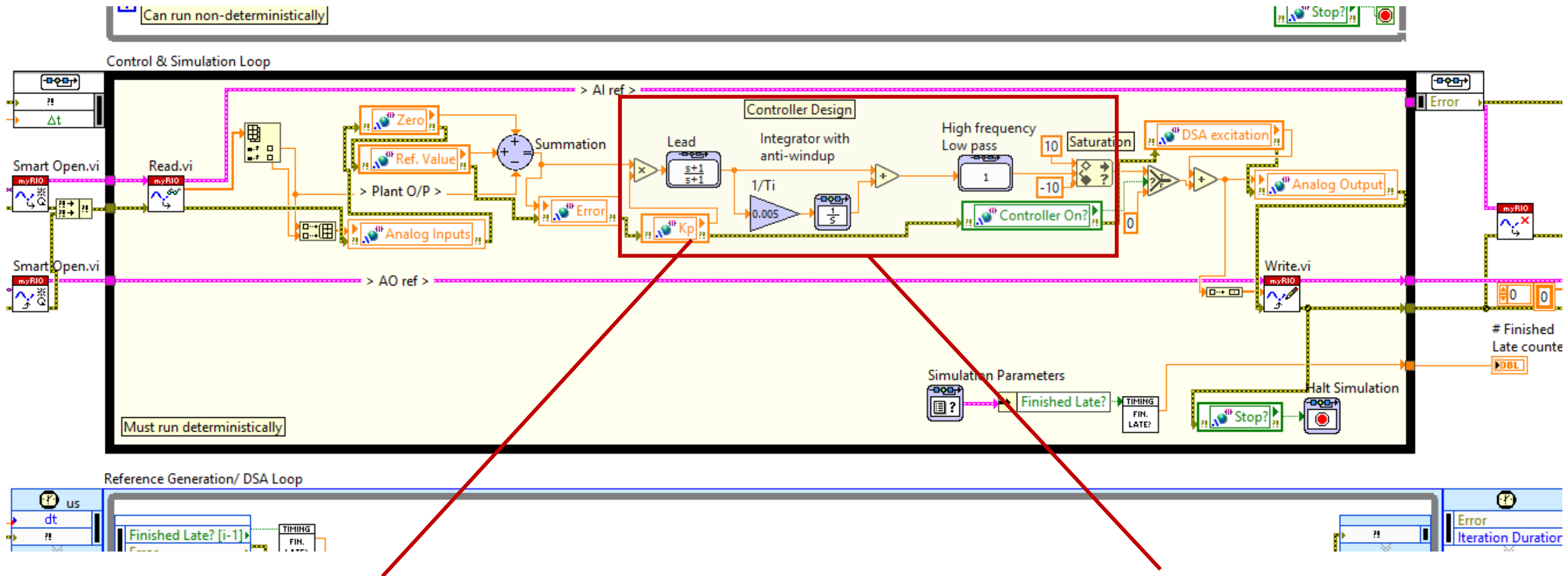

Proportional gain value comes from the control on the block diagram.

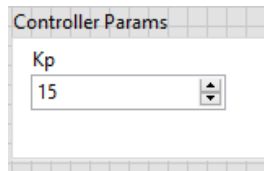

Controller template with Proportional, Lead, integrator with anti-wind up and high frequency roll-off.

Double-click on each block to change its value. Press Ctrl+H for in-context help.


Using the Dynamic Signal Analyzer

1. Select options on the User interface:

- Select DSA tab in the Waveform control.
- Select DSA tab in the Ref. Gen. control.

Select frequency response:

- **Frequency response**
 - **Plant** – Plant, takes AO 0 (control effort) as its input and gives AI 0 as its output (measured plant output).
 - **RR (Loop Return Ratio)** –
Input: Error
Output: Plant output
 - **CL (Closed loop)** –
Input: Reference
Output: Plant output

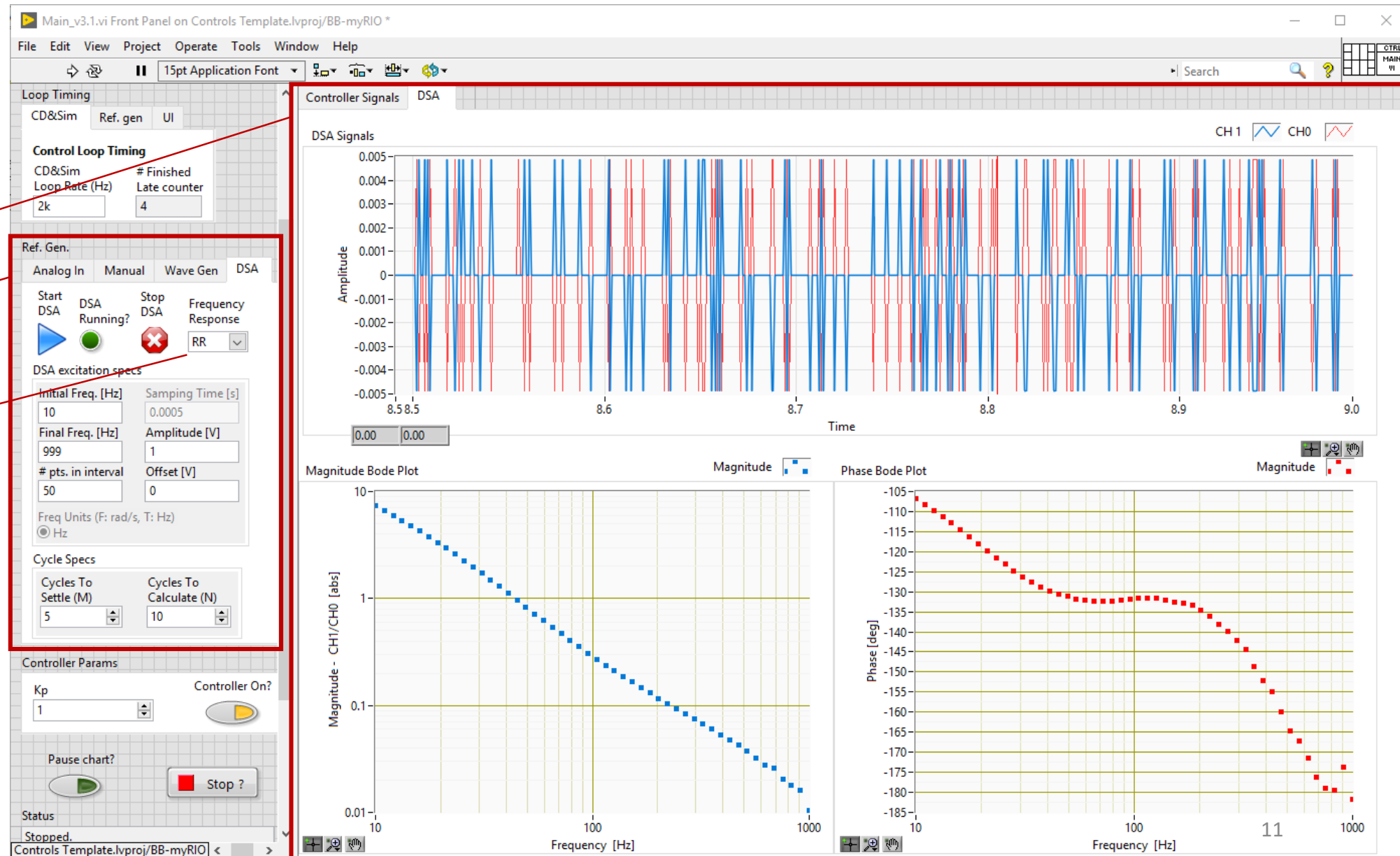


Using the Dynamic Signal Analyzer

2. Select DSA options:

- **Initial frequency** – low frequency can take a long time to capture.
- **Final frequency** – maximum limited to Nyquist frequency of the loop sample time ($f_s/2$)
- **# pts. in interval** – number of points logarithmically distributed between initial and final frequency.
- **Amplitude** – input signal at which to get the frequency response.
- **Offset** – to add a bias to the input signal.
- **Cycles to settle (M)** – Number of sine cycles for setting at each frequency before measurement starts.
- **Cycles to calculate (N)** – Number of sine cycles to average over to make measurement at a frequency

Note: For low frequency, use low values of M and N to reduce measurement time (2,2). For high frequency use large M and N to get more stable measurements (20,50).

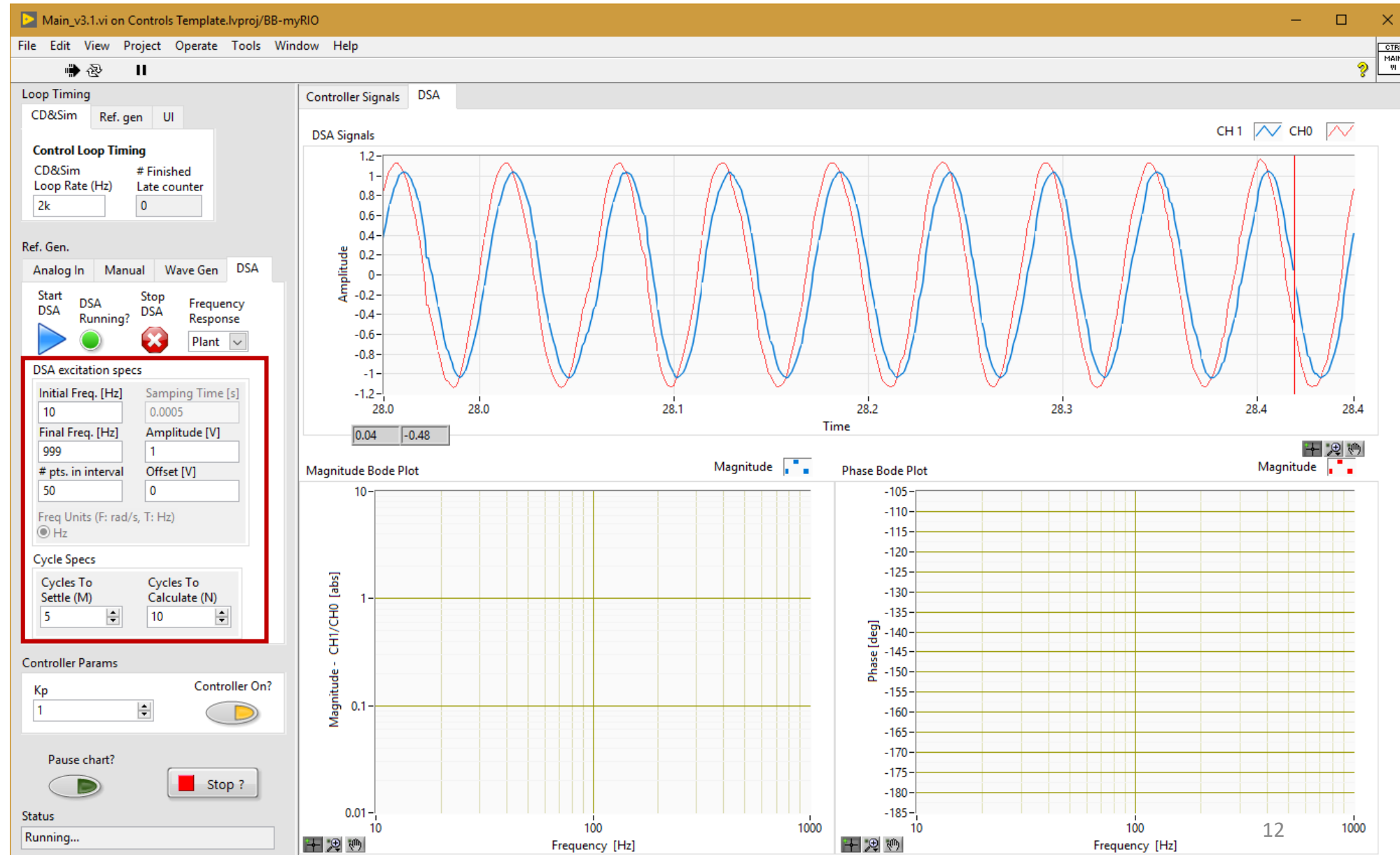


Using the Dynamic Signal Analyzer

3. Start the DSA:

- Press start to start the DSA.
- When running, the DSA running indicator will be on.
- It can take a while, especially if starting from low frequencies.
- Press stop to abort the Capture.


Note:

- Make sure that your control effort doesn't saturate as you sweep the frequencies. If so, reduce amplitude.
- If amplitude of error signal for the RR node plot is too small, reduce Kp, get the RR node plot, export it to excel and then manually shift the magnitude plot (abs) upwards by Kp times.

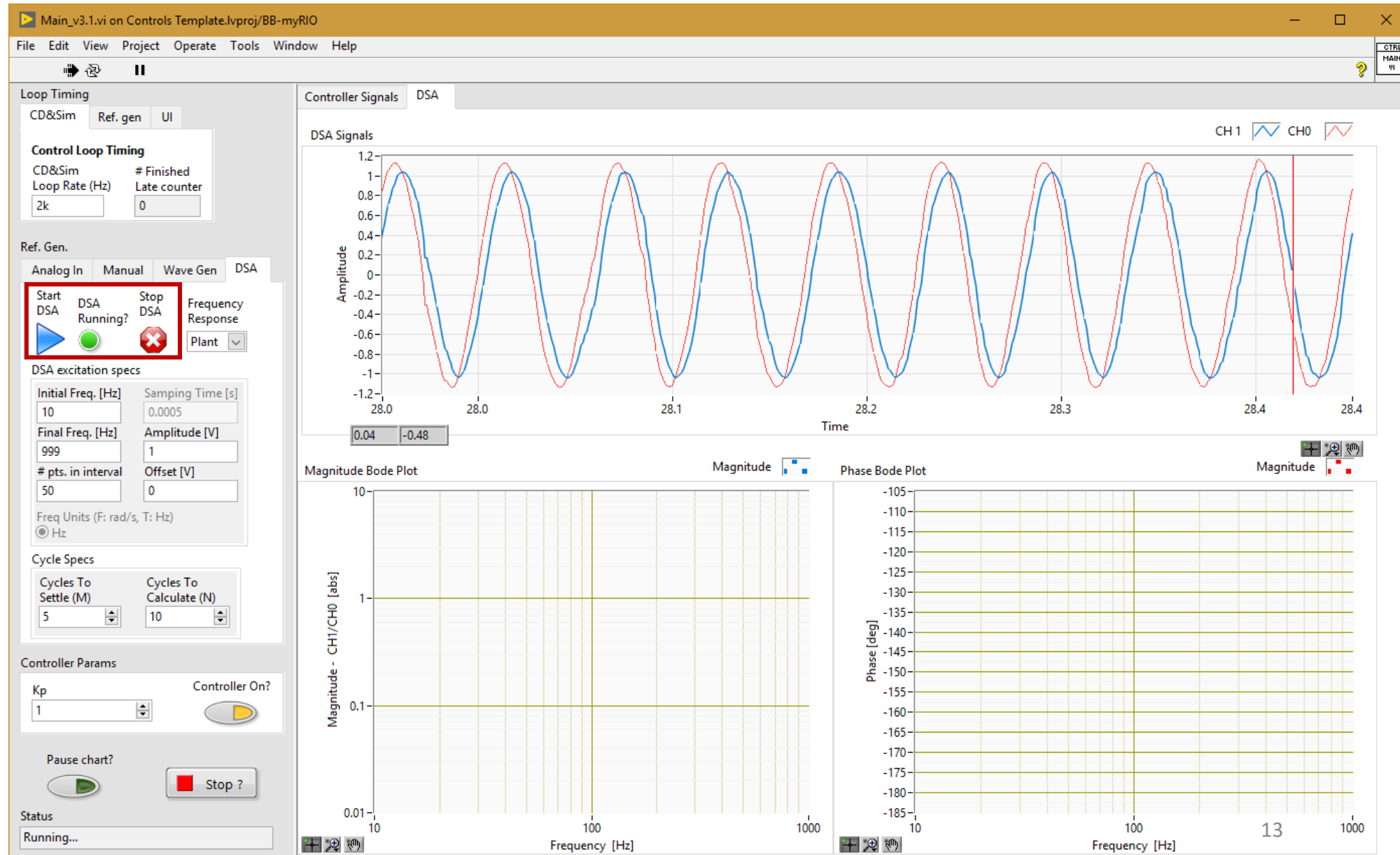


Using the Dynamic Signal Analyzer

4. **Export the Mag and phase bode plots**
- Once the DSA stops, it will plot the magnitude and phase plots.
 - You can export the magnitude and phase plots to excel by right-click and selecting Export->to Excel as shown.

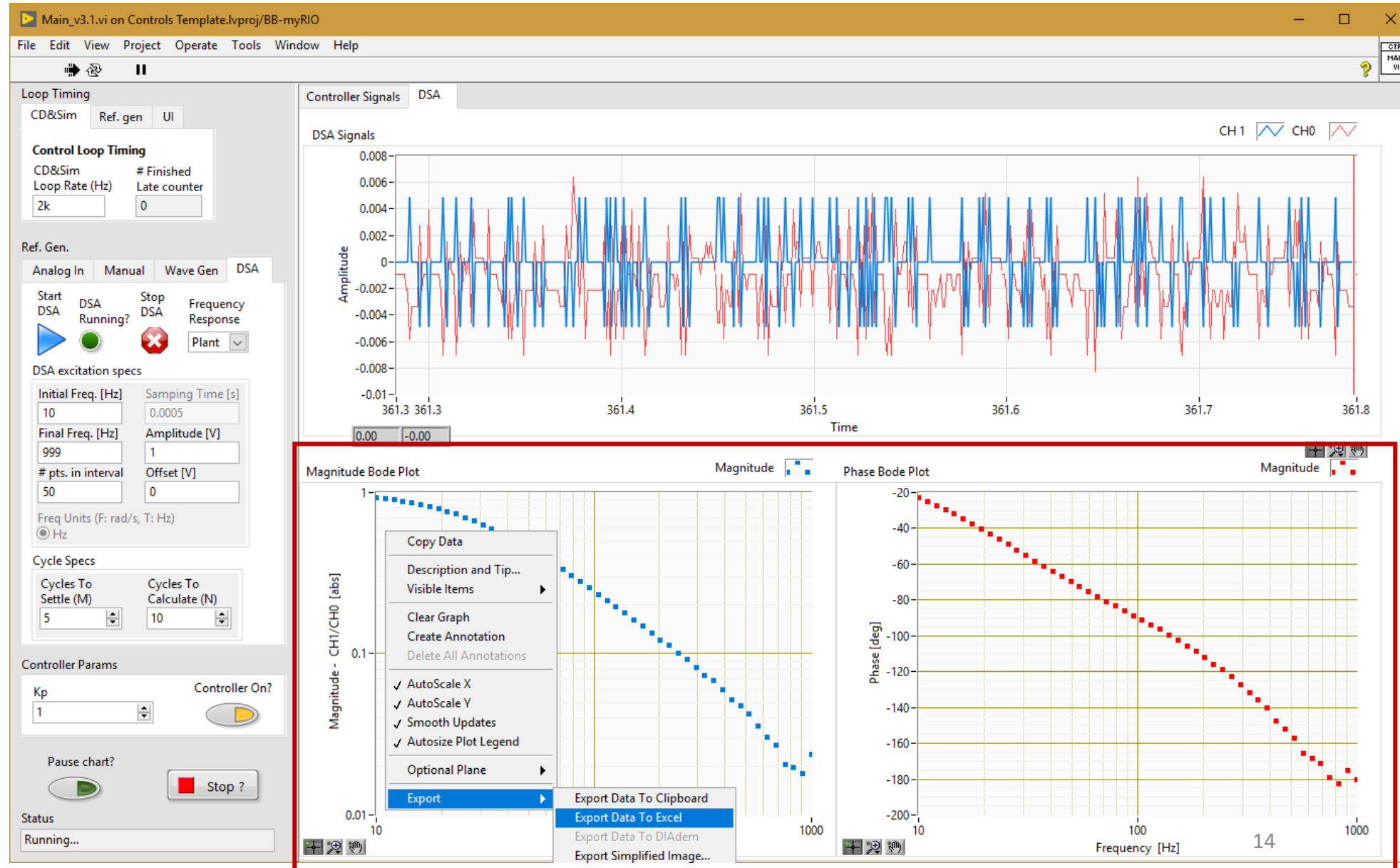

Sample DSA results

Unity plant (AO connected to AI), $T_s = 0.5$ msec

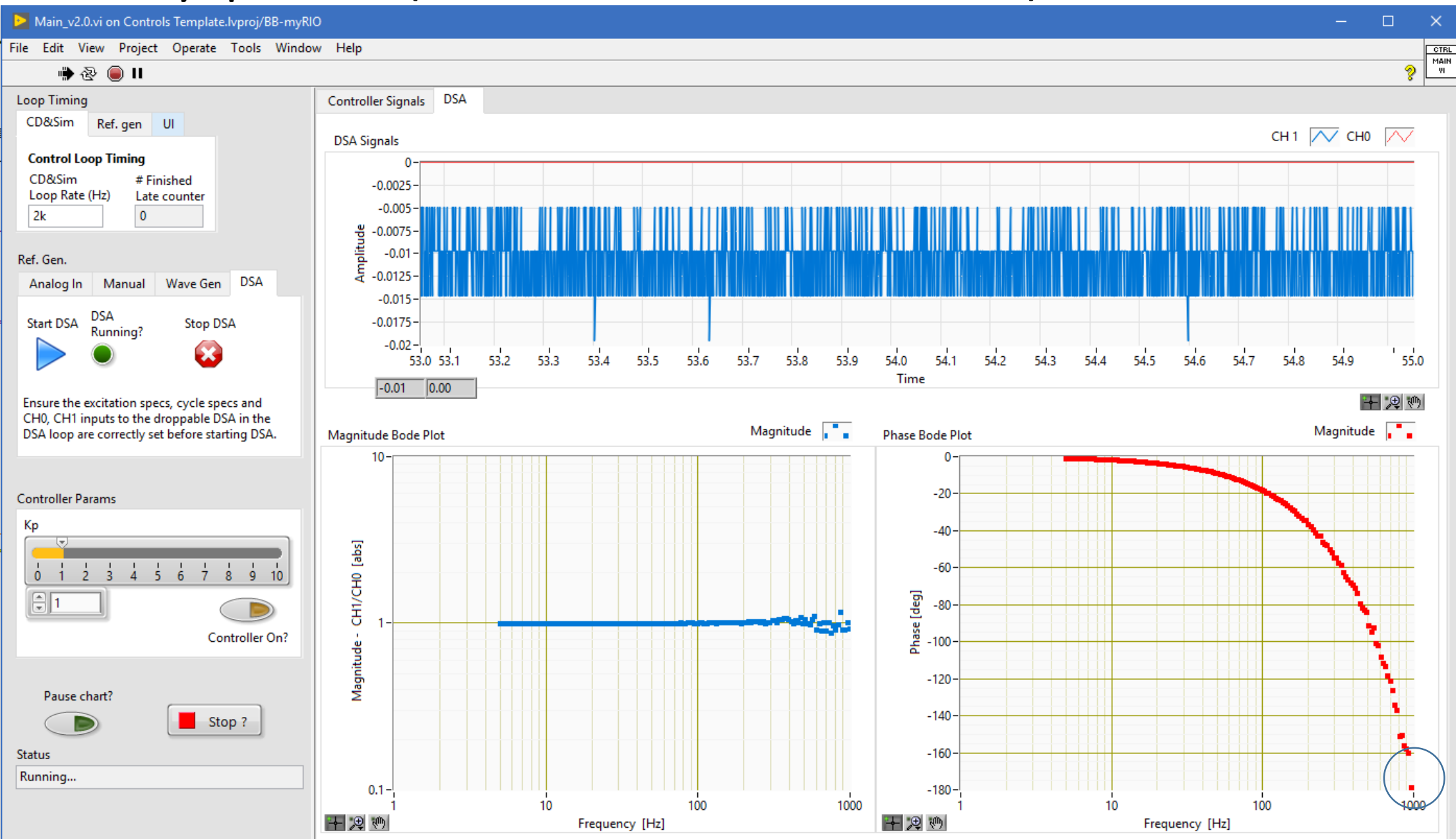


Transfer function in continuous time
 $H(s) = 1$

As seen by the myRIO, due to A/D sampling, there is a one sample delay,
 $H(s) = 1 \cdot e^{-s \cdot T_s}$
i.e.,
 $H(z) = z^{-1}$


Thus, one loop sample time delay.

$f_s = 2$ kHz
180 deg. Phase lag at Nyquist freq. (1 kHz)

RC low pass filter, $\tau = 4.7 \text{ ms}$, $f_{3dB} \approx 34 \text{ Hz}$

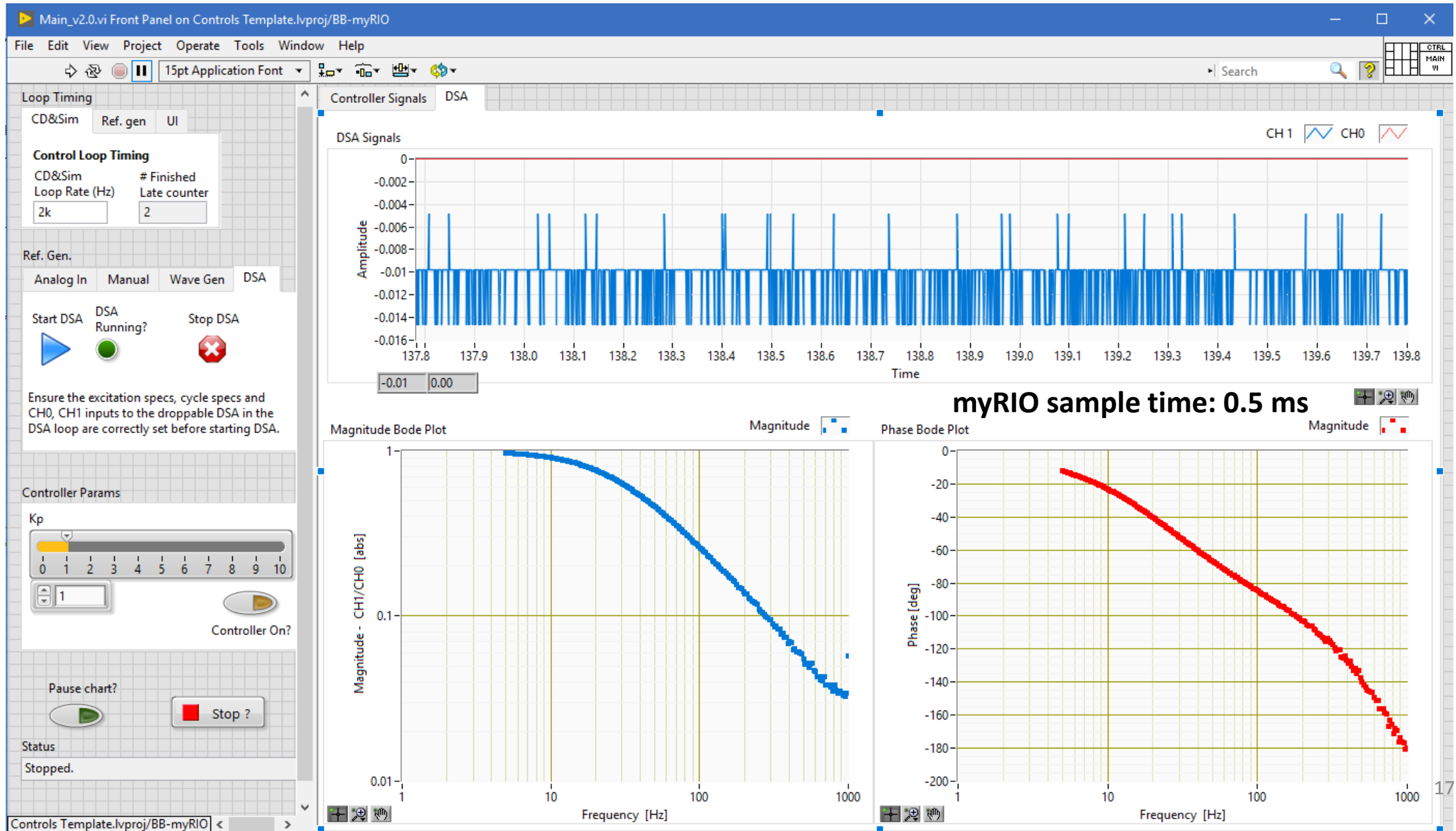

Another low pass filter, $\tau = 6.2 \text{ ms}$, $f_{3dB} \approx 25 \text{ Hz}$, $T_s = 0.0005 \text{ s}$

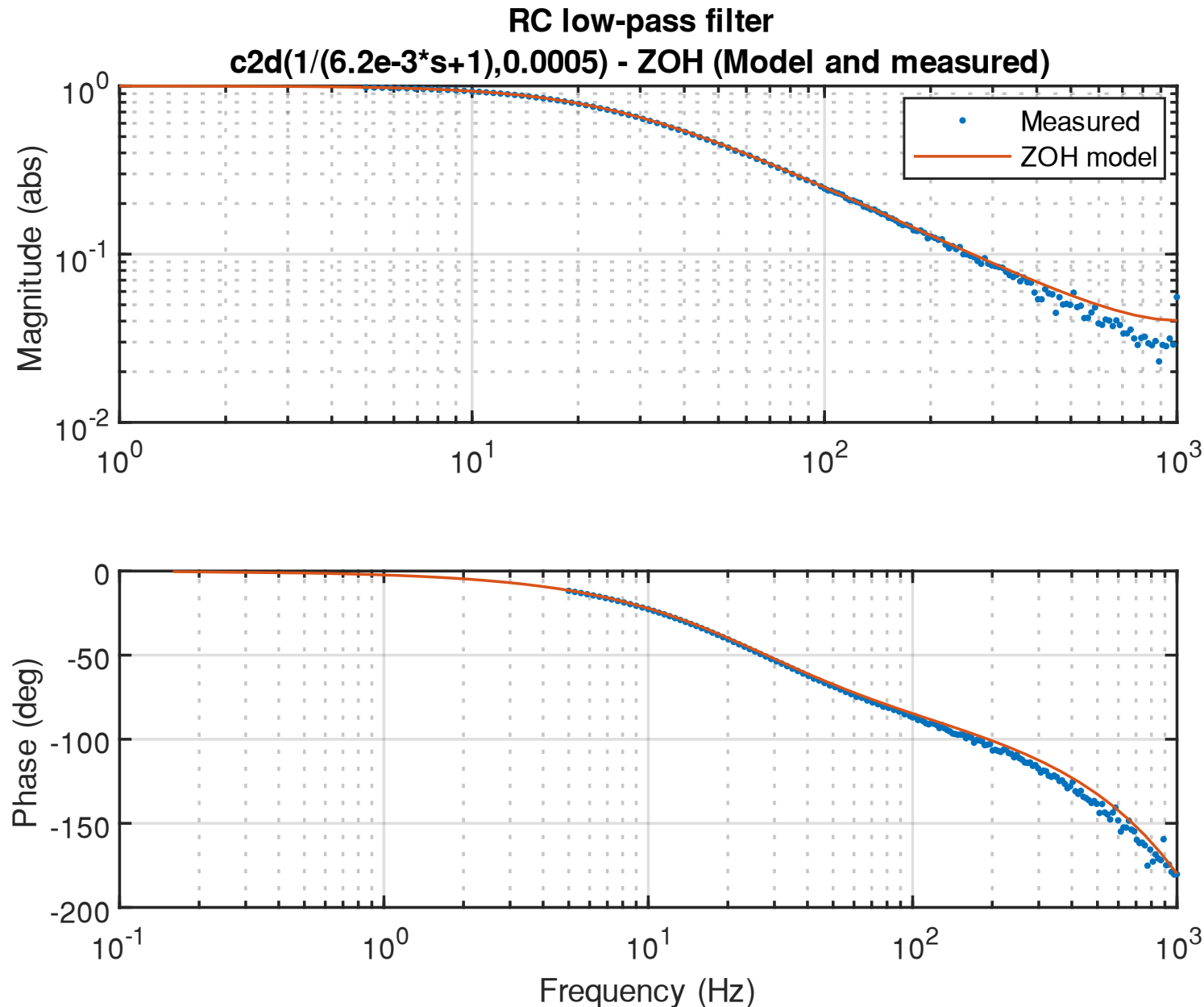

In MATLAB, to convert from continuous time transfer function $H(s)$ to discrete time transfer function $H(z)$, use the `c2d` command:


$$\text{sys_d} = \text{c2d}(\text{sys_c}, T_s)$$

T_s = sample time, sys_c = continuous time transfer function and sys_d = discrete time transfer function.

We obtain the zero-order hold equivalent of our transfer function in discrete domain. That will match pretty well with the bode plots you obtain from the DSA functionality of the program.

More information

Excitation specs and cycle specs

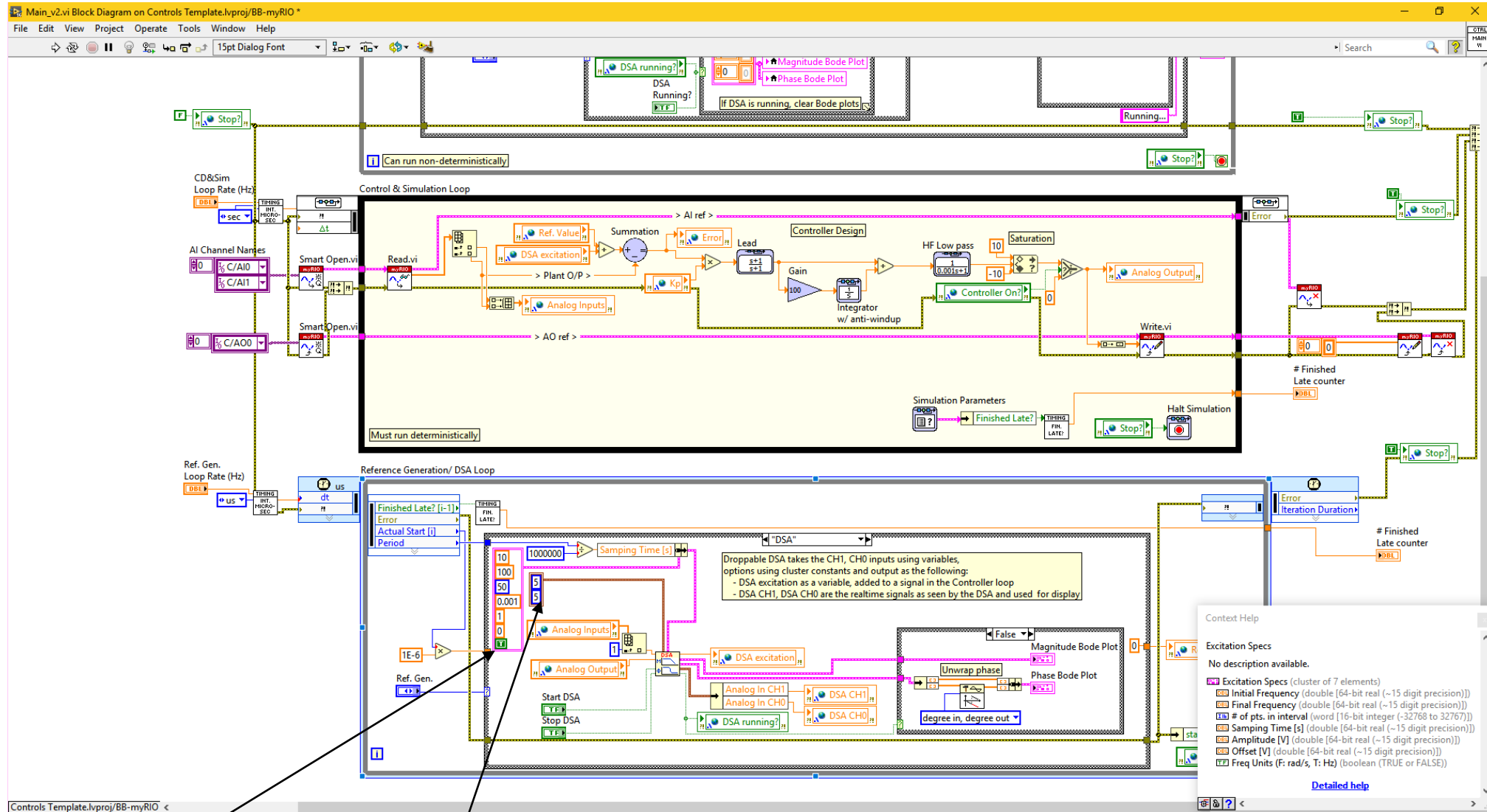


Cycle specs are directly entered as constants into the Droppable DSA. Hover over them with context help on (Ctrl+H) to see the description.

Reference value set to 0 when the DSA mode is selected.




DSA running

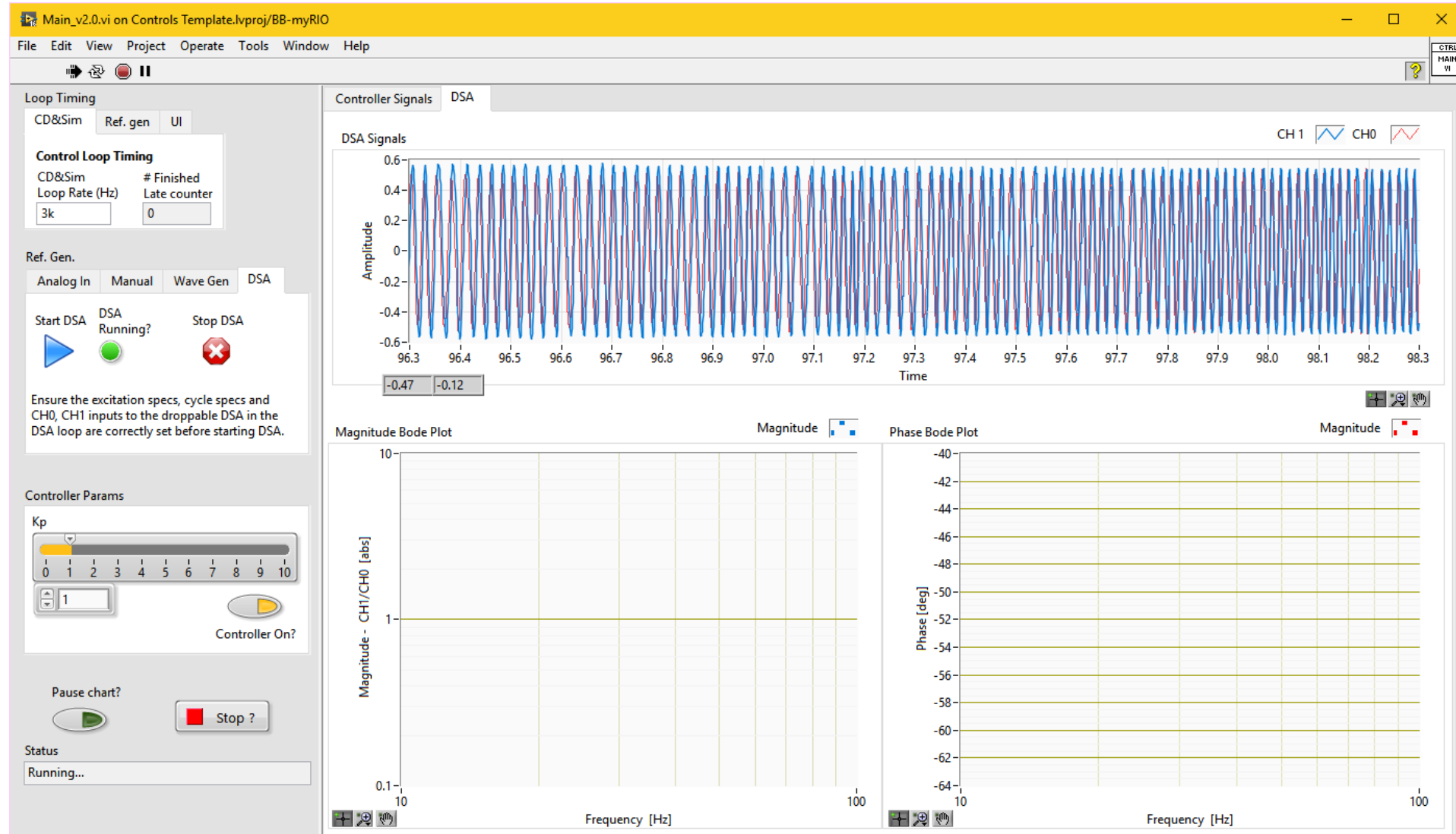

When DSA is running the Bode plot is cleared and it is replotted with the DSA generated data after it stops.

DSA results

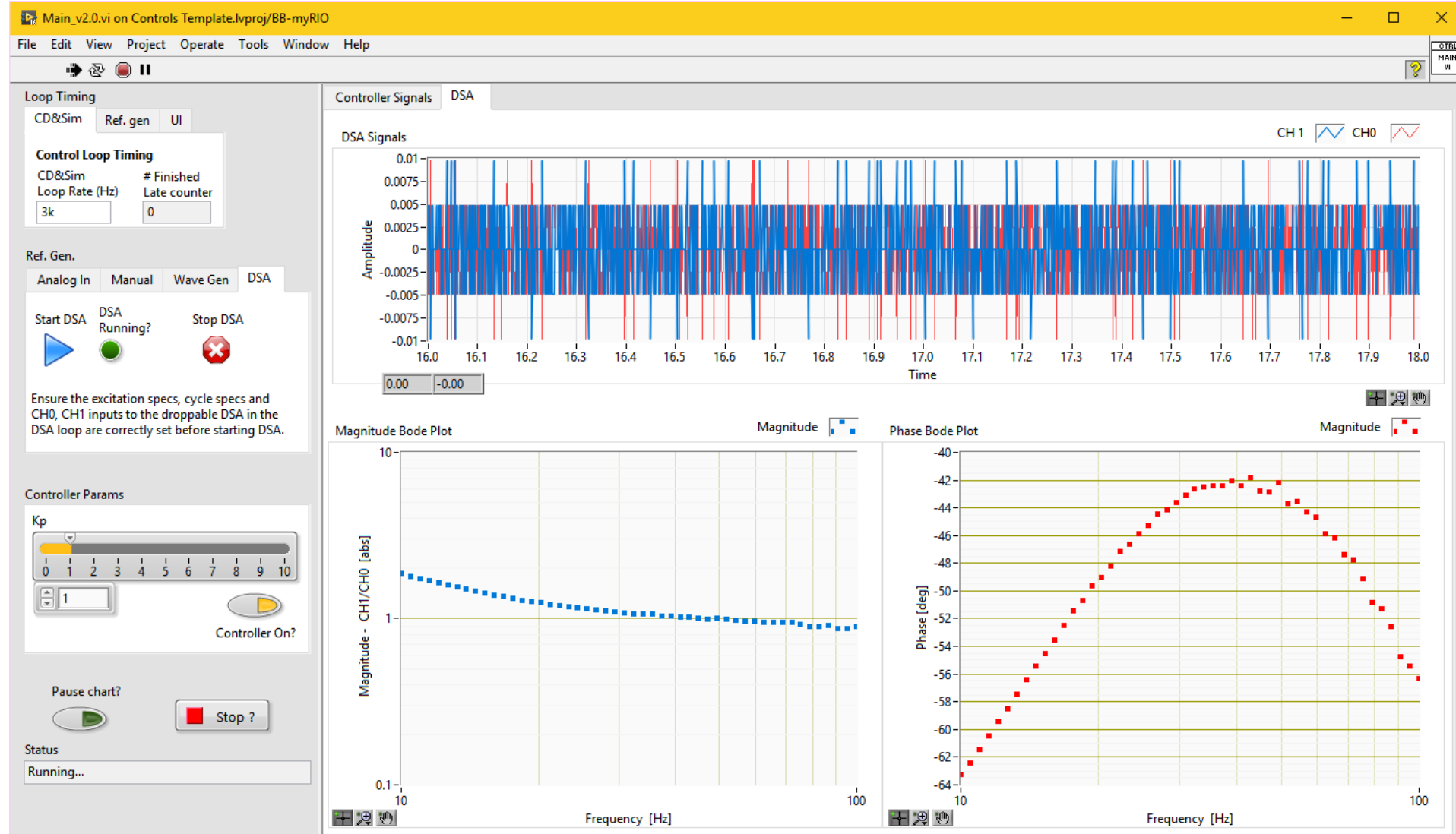

Note: There is a time delay of one sample time (T_s) due to the DAQmx implementation in which it sends the Analog Output for this sample at the clock edge for the next sample. (<http://www.ni.com/tutorial/3215/en/>, <http://zone.ni.com/reference/en-XX/help/370466Y-01/mxcncpts/controlappcase6/>)